

enbio®

Enbio PRO. Small footprint. Big space.



Triple Vacuum
Full Drying

Designed  in Switzerland



Triple Vacuum
Full Drying

Enbio PRO

S / 273°F
9 min

Unwrapped

B / 273°F
15 min


Wrapped
(pouches)

B / 250°F
45 min

Delicate
instruments



16 diagnostic sets OR
16 handpieces OR
1 – 4 instrument cassettes OR
1 implant cassette

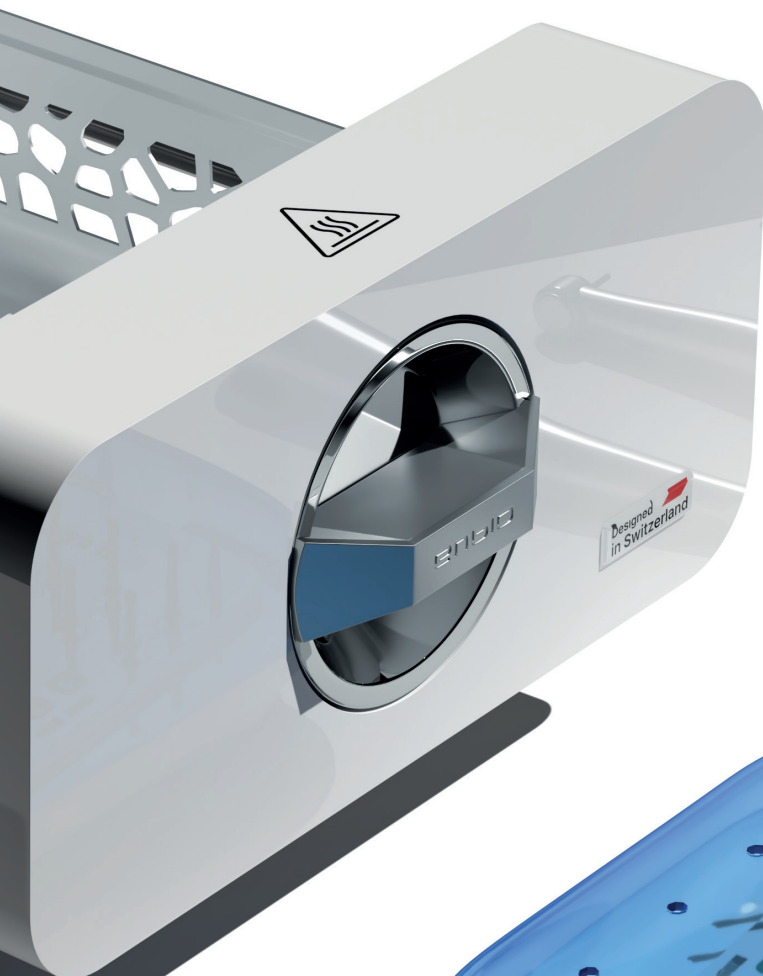

Enbio PRO. Small footprint. Big space.

Workhorse in a small practice.
Double capacity and space
for implant and surgical cassettes.



Usable load space
11.2 × 7.5 × 2.8 in

16 diagnostic sets OR
16 handpieces OR
1–4 instrument cassettes OR
1 implant cassette

small instrument tray
pouch rack



Enbio S. Speedmaster.

Unbeatable speed in the class.
15 minute vacuum cycle.
Bone dry pouches.



Usable load space

11.2 × 7.1 × 1.3 in


8 diagnostic sets OR
8 handpieces




Performance technology inside.

Die-cast Aluminum Chamber

Patented aluminum alloy chamber is manufactured through an advanced die-casting process that yields exceptional structural integrity and thermal conductivity, ensuring uncompromised performance throughout the device's operational lifespan.

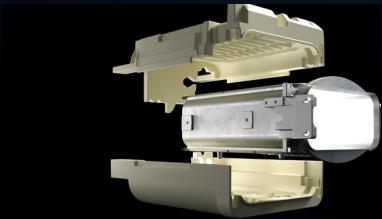

**Polyurethane foam
chamber insulation**

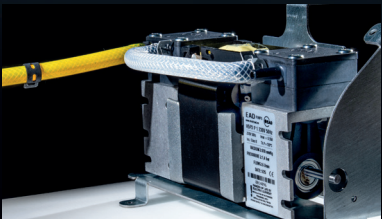

**Heavy-duty
vacuum pump**



Lifetime warranty

The chamber is covered by a lifetime warranty.

Structure-integrated heaters

Heating elements are directly integrated into the die-cast chamber structure during production. Seamless thermal interface significantly reduces warm-up intervals while ensuring optimal temperature uniformity.




Fully automatic with touch screen control



Red silicone chamber seal with snap-off tray mount

Enbio PRO

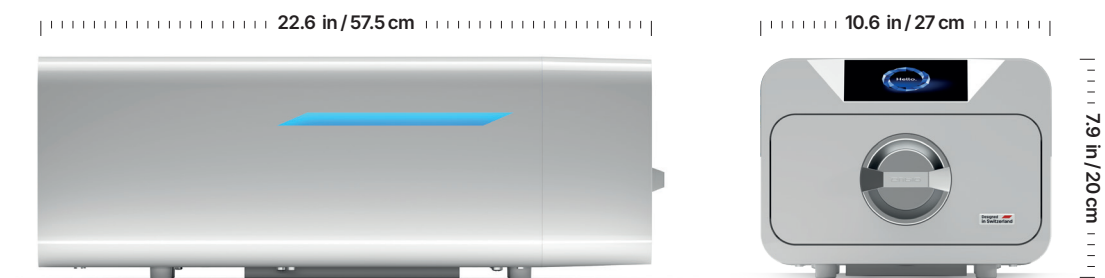

Usable load space
11.2 × 7.5 × 2.8 in

Unit weight

44 lbs / 20 kg

Capacity

Max. packages amount (circa)	16
Max. recommended batch weight	3.5 lbs / 1.6 kg
Process chamber volume	5.3 l

Electrical

Power supply	110-120V AC/60 Hz
Power	2.45 kW
Max. power consumption	15 A

Pressure & temperature

Working pressure	2.1 Bar / 30.47 psi
Max. pressure	2.3 Bar / 33.38 psi
Max. process temperature	137°C / 278.6°F

Misc

Ingress protection	IP20
Data archiving	USB stick
Water conductivity*	< 15 µS/cm

Medical device class II b
EN-13060
93/42/EEC Directive

FDA*
Cleared

TUVNORD

GXE CLEANING
STERILIZATION

MDSAP
Medical Device
Single Audit Program



Enbio Corp.
1075 HWY 124, Hoschton G.A. 30548, Bldg 300 Ste 300

enbio.com