

 symmetry surgical® | **Bovie**®  
**HIGH FREQUENCY DESICCATORS**



# BOVIE® DERM 942 HIGH FREQUENCY DESICCATOR




The DERM 942 and its patent-protected handpiece are designed and engineered to be the most reliable and durable desiccator available today. The unit's life is prolonged by incorporating the latest in chip technology and reducing the number of internal parts, thereby generating less internal heat.

The DERM 942 comes with 40 watts of coagulation capabilities and is adjustable in 1/10th watt increments in the first 10 watts (in monopolar and bipolar). It is a ground reference unit which means no grounding pad is required. The DERM 942 is the most cost-effective and user-friendly unit on the market, meaning less to do and less to buy. Set-up and usage are very simple and straightforward when compared to the competitor. For example, adjusting from low to high power and adding foot control to the handpiece are accomplished in the following ways.

Bovie DERM 942	Competitor
<b>Adjusting from Lo to Hi Power</b>	
<b>1-Step</b> <b>Adjust</b> rotary dial or press handpiece buttons	<b>4-Steps</b> <b>Unplug</b> handpiece from Lo power receptacle <b>Move</b> power selector switch from Lo to Hi <b>Plug</b> handpiece into Hi power receptacle <b>Adjust</b> rotary dial or press handpiece buttons
<b>Adding Foot Control to Handpiece</b>	
<b>1-Step</b> <b>Plug</b> footswitch into front of unit	<b>6-Steps</b> <b>Unplug</b> handpiece from bottom of unit <b>Unplug</b> handpiece from front of unit <b>Detach</b> handpiece from holder <b>Attach</b> foot control handpiece to holder <b>Plug</b> foot control handpiece into front of unit <b>Plug</b> footswitch into bottom of unit

# BOVIE® DERM 942 HIGH FREQUENCY DESICCATOR

## PRODUCTS INCLUDED

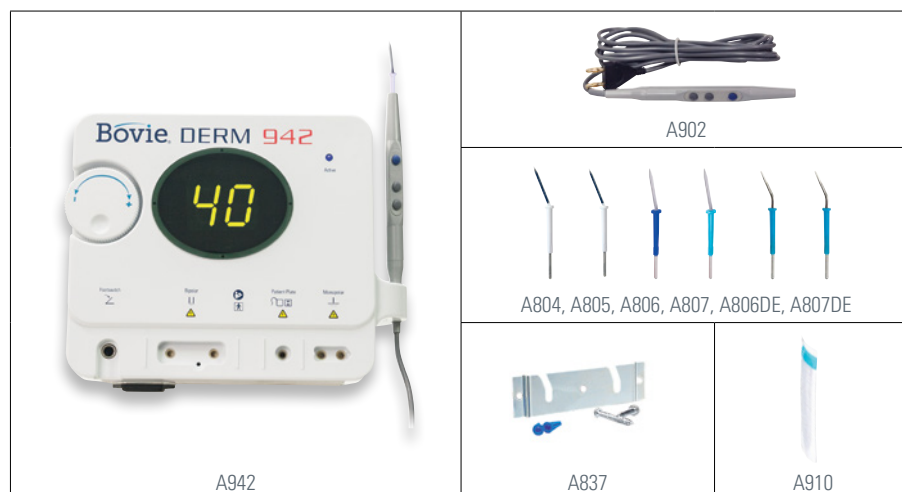


## Style # A942

### ▼ Products Included

- A942** (1) high frequency desiccator
- A902** (1) handpiece
- A804** (5) non-sterile sharp dermal tips
- A805** (2) sterile sharp dermal tips
- A806** (5) non-sterile blunt dermal tips
- A807** (2) sterile blunt dermal tips
- A806DE** (2) non-sterile DERM-Elite™ Premium blunt dermal tips
- A807DE** (2) sterile DERM-Elite Premium blunt dermal tips
- A910** (2) non-sterile disposable handpiece sheaths
- A837** (1) wall mount kit, 110 VAC Hospital-grade power cord\*  
User's guide (**MC-55-239-001**)\*

## ▼ COMPLEMENTARY PRODUCTS

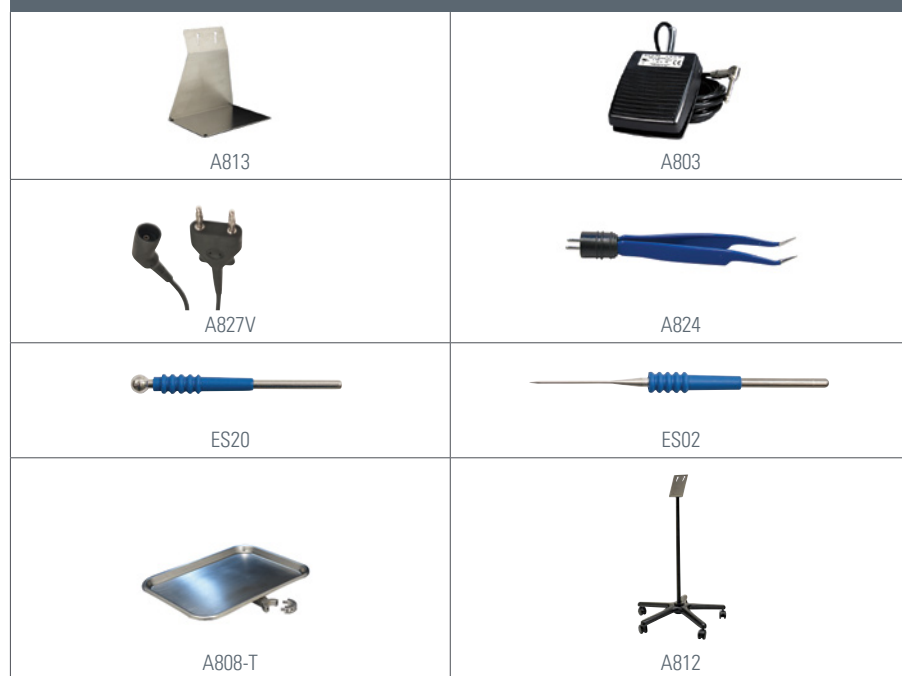

### Dimensions

Height: 188mm (7.40")  
Width: 228mm (8.98")  
Depth: 105mm (4.13")  
Weight: < 2.26kg (< 5lbs.)

### Output Power

To 40 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

### Line Voltage

100-240 VAC

**Output Frequency:** 328-492 kHz

**Line Frequency:** 50-60 Hz

**Warranty:** FOUR YEARS

\*Not Pictured

†Sold separately

# BOVIE® BANTAM | PRO ELECTROSURGICAL GENERATOR



The Bovie® Bantam PRO A952 is a high frequency desiccator that also has cutting ability. The unit provides 50 watts of cutting power and coagulation modes – plenty of power for most procedures performed in the doctor's office. Power is adjustable in 1/10th watt increments from 0 to 10 watts, for the most accurate output available. No grounding pad is required when working in desiccation mode.

The unit also features NEM™ (Neutral Electrode Monitoring) pad sensing technology which monitors the return electrode for optimum patient protection. The Bovie® Bantam PRO A952 is backed by Bovie's four (4) year manufacturer's warranty.

# BOVIE® BANTAM | PRO ELECTROSURGICAL GENERATOR

## PRODUCTS INCLUDED

 <p>A952</p>	
 <p>A902</p>	 <p>A804, A805, A806, A807, A806DE, A807DE</p>
 <p>ES01, ES02, ES20</p>	 <p>A837</p>
 <p>A910</p>	 <p>A910ST</p>

## Style # A952

### ▼ Products Included

- A952** (1) electrosurgical generator
  - A902** (1) handpiece
  - A804** (5) non-sterile sharp dermal tips
  - A805** (2) sterile sharp dermal tips
  - A806** (5) non-sterile blunt dermal tips
  - A807** (2) sterile blunt dermal tips
  - A806DE** (2) non-sterile DERM-Elite™ Premium blunt dermal tips
  - A807DE** (2) sterile DERM-Elite Premium blunt dermal tips
  - ES01** (3) blade sterile disposable electrodes
  - ES02** (1) needle sterile disposable electrodes
  - ES20** (1) ball sterile disposable electrodes
  - A910** (2) non-sterile disposable handpiece sheaths
  - A910ST** (2) sterile disposable handpiece sheaths
  - A1252C** (2) reusable grounding cord\*
  - ESRE** (5) disposable split grounding pads\*
  - A837** (1) wall mount kit
- 110 VAC Hospital-grade power cord\*  
User's guide (MC-55-238-001)\*

### Dimensions

Height: 188 mm (7.40 in.)  
Width: 228 mm (8.98 in.)  
Depth: 105 mm (4.13 in.)  
Weight: < 2.26 kg (< 5 lbs.)

### Output Power

Cut: 50 watts  
Blend: 50 watts  
Coagulation: 50 watts  
Fulguration: 40 watts  
Bipolar: 40 watts  
Micro Bipolar: 40 watts

**Output Frequency:** 304-451 kHz

**Line Frequency:** 50-60 Hz

**Line Voltage:** 100-240 VAC

**Warranty:** **FOUR YEARS**









# BOVIE® BANTAM | PRO A952-G ELECTROSURGICAL GENERATOR



The Bovie® Bantam PRO | A952-G is packaged with the essential accessories needed to perform an OB/GYN LEEP procedure. The full system is situated on a mobile stand for effortless portability. It is an extremely affordable electrosurgical generator package with smoke evacuation and includes a variety of electrodes and LLETZ loop electrodes. The unit is backed by Bovie's four (4) year manufacturer's warranty.

# BOVIE® BANTAM | PRO A952-G ELECTROSURGICAL GENERATOR

## PRODUCTS INCLUDED

 <p>A952</p>			
 <p>A902</p>	 <p>A1252C</p>	 <p>SE02</p>	
 <p>A804, A805, A806, A807, A806DE, A807DE</p>		 <p>A910ST</p>	 <p>A812C</p>
 <p>ES01, ES02, ES20</p>	 <p>ESLK</p>	 <p>ES07</p>	
 <p>786T</p>	 <p>A837</p>	 <p>ESRE</p>	 <p>A910</p>
 <p>SERF</p>			

## Style # A952-G

### ▼ Products Included

#### ▼ Electrosurgery Equipment

- A952** (1) electrosurgical generator
  - A902** (1) handpiece
  - A804** (5) non-sterile sharp dermal tips
  - A805** (2) sterile sharp dermal tips
  - A806** (5) non-sterile blunt dermal tips
  - A807** (2) sterile blunt dermal tips
  - A806DE** (2) non-sterile DERM-Elite™ Premium blunt dermal tips
  - A807DE** (2) sterile DERM-Elite™ Premium blunt dermal tips
  - ES01** (3) blade sterile disposable electrodes
  - ES02** (1) needle sterile disposable electrodes
  - ES20** (1) ball sterile disposable electrodes
  - A910ST** (2) sterile disposable handpiece sheaths
  - A910** (2) non-sterile disposable handpiece sheaths
  - A1252C** (2) reusable grounding cord
  - ESRE** (5) disposable split grounding pads
  - A837** (1) wall mount kit
  - ES07** (5) extended ball electrodes
- 110 VAC Hospital-grade power cord\*  
User's guide (**MC-55-238-001**)\*

#### ESLK electrosurgical loop kit

- ES09** (1) 10mm x 10mm loop
- ES11** (1) 15mm x 10mm loop
- ES12** (1) 20mm x 8mm loop
- ES13** (1) 20mm x 15mm loop
- ES16** (1) 10mm x 8mm square

#### A812-C Mobile stand

- A812** (1) mobile stand for A952
- A812-E** (1) height extender for A812
- A808-T** (1) top tray and clamp for A812
- A812-BT** (1) bottom tray for A812

#### Smoke evacuation system

- SE02** (1) Smoke Shark® II with pneumatic footswitch
  - SF35** (1) 35-hour filter
  - 786T** (1) 7/8" (22.2mm) tube x 6' (182.9cm) long
  - SERF** (1) reducer fitting (attaches to vaginal speculum)
- 110 VAC Hospital-grade power cord\*  
User's guide (**MC-23088**)\*

# BOVIE® DERM 102 HIGH FREQUENCY DESICCATOR



The Bovie® DERM 102 High Frequency Desiccator offers ten watts of monopolar power, adjustable in 1/10th watt increments, to desiccate and fulgurate tissue. The power can be adjusted using the 3-Button Power Control Handpiece A902 (Patent #8,998,899), or using the rotary dial on the unit. The unit features a large, bright, LED display for excellent visibility across the procedure room. If foot control activation is preferred while using the DERM 102, just plug in the A803 Bovie Footswitch (sold separately) and activate. No separate handpiece is required.

For the practitioner who wants to also use bi-polar electrosurgery along with high frequency desiccation, the DERM 102 High Frequency Desiccator provides 10 watts of bipolar power, adjustable in 1/10th watt increments. To utilize this feature, the A827V\* bipolar cable, a bipolar forcep\* and footswitch activation are required.

The DERM 102 is packaged ready to use. The unit comes with a total of 18 disposable dermal tips. By the time the practitioners utilizes all 18 electrodes, they will have paid for the unit (twice), and generated income for their practice.



# BOVIE® DERM 102 HIGH FREQUENCY DESICCATOR

## PRODUCTS INCLUDED



## Style # DERM 102

### ▼ Products Included

- DERM 102** (1) high frequency desiccator  
**A902** (1) handpiece  
**A910** (2) non-sterile, disposable handpiece sheaths  
**A804** (5) non-sterile single-use sharp dermal tips  
**A805** (2) sterile single-use sharp dermal tips  
**A806** (5) non-sterile single-use blunt dermal tips  
**A807** (2) sterile single-use blunt dermal tips  
**A806DE** (2) non-sterile DERM-Elite™ Premium blunt dermal tips  
**A807DE** (2) sterile DERM-Elite Premium blunt dermal tips  
**A837** (1) wall mount kit  
 110 VAC Hospital-grade power cord\*  
 User's guide CD (**MC-55-229-001**)\*

### Dimensions

Height: 208mm (8.2")  
 Width: 184mm (7.2")  
 Depth: 69mm (2.7")  
 Weight: < 1.4kg (< 3 lbs.)

### Output Power

Up to 10 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

### Line Voltage

100-240 VAC

### Line Frequency: 50 – 60 Hz

**Main Current:** 320 mA Max

### Warranty: **FOUR YEARS**

\* Not Pictured



3034 Owen Drive | Antioch, TN 37013 | 800.251.3000 (US) | +41 (0) 52 557 59 00  
customerservice@symmetrysurgical.com | [symmetrysurgical.com](https://www.symmetrysurgical.com)  
©2022 Symmetry Surgical Inc. All rights reserved.  
Bovie High Frequency Dessicator Brochure 03.23.2022 v1