


## PowerSpin™ MH Centrifuge

### Product Images

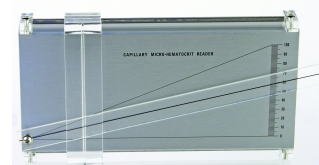

### Short Description

POWERSPIN™ MH MICRO-HEMATOCRIT CENTRIFUGE  
FIXED SPEED 12000 RPM  
24 PLACE MICROHEMATOCRIT ROTOR  
5 MIN MECHANICAL TIMER  
POWER INPUT 110V  
MICROHEMATOCRIT READER INCLUDED.

## Description

Model C-MH30 is made of ABS plastic with inner metal chamber. A safety switch built into the handle assembly disconnects the power to the motor whenever the latch is lifted.

The 24 place microhematocrit rotor is an aluminum shell for durability and lighter weight

The brushed motor is mounted on rubber mounts to provide quieter running and less vibration. cleaned and disinfected.

## Additional Information

Accessories	Header	Accessories
	CMH30-Reader	Microhematocrit reader
	CMH30-03	Gasket Replacement For CMH30 Centrifuge
Specification	Rotor	24 place
	Power Input	110V
Size and Weight	Size	11.6 in x 10.3 in x 9.3 in (295 x 261 x 235 mm)
	Shipping Size	12.7 in x 12.7 in x 12.2 in (323 x 323 x 311 mm)
	Weight	17.4 lbs (7.9 kg)
	Shipping Weight	19.4 lbs (8.8 kg)